

液相色谱法分析脂溶性维生素 Analysis of Fat-soluble Vitamin

► 分析条件

Column : ZORBAX XDB-C8(4.6*150mm)


UV/Vis Wavelength : 280nm

Mobile phase : A=Water B=MeOH (A:B=5:95)

Flow rate : 1ml/min

Inject volume: 10ul

► 色谱图

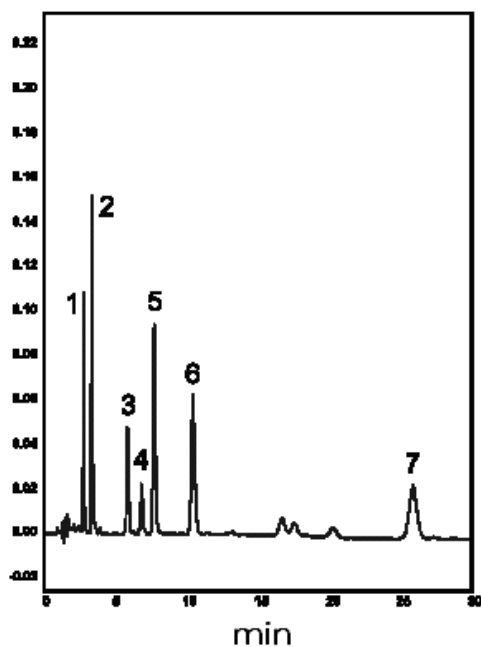

1. Retinol/维生素A
2. Retinol Acetate/醋酸维生素A
3. Vitamin D3/维生素D3
4. α -Tocopherol/ α -维生素E
5. β -Tocopherol/ β -维生素E
6. Tocopherol Acetate/醋酸维生素E
7. Retinol Palmitate



Younglin 氨基酸分析系统